


# Installation Schematic for *VistaClear*™ MODEL 1000

Pat. No. 6,423,219 and other patents pending

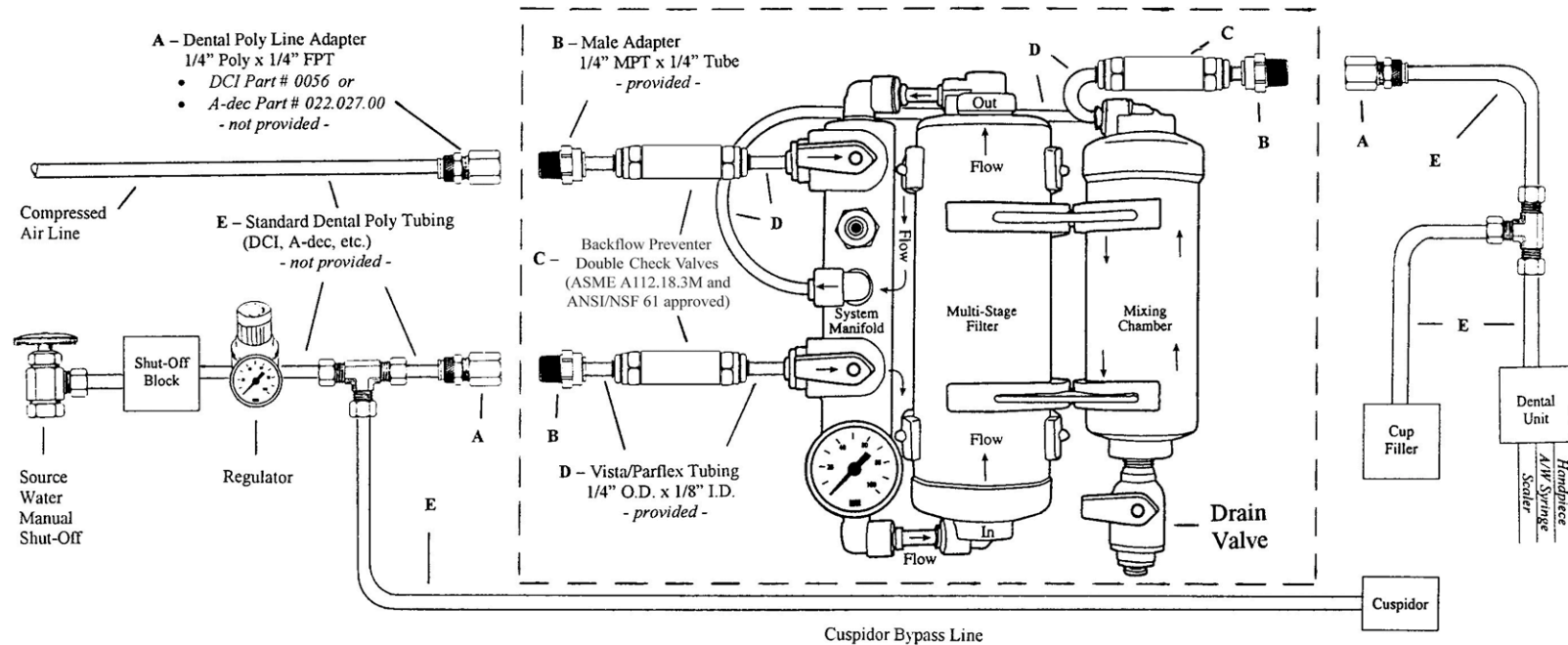

## IMPORTANT

- All items within the dashed lines are included with the *VistaClear* system (except mounting screws)
- *Vista (Parflex)* tubing **MUST** be used on the system and all fittings within the dashed lines (1/4" O.D. x 1/8" I.D. tubing is included)
- Do **NOT** use tubing inserts with *Vista (Parflex)* tubing
- The *VistaClear* system **MUST** be connected to the water system **AFTER** the dental waterline shut-off block and pressure regulator
- If a cuspidor is present, **ALWAYS** feed it directly from the unfiltered water supply after the regulator – **NEVER** from the *VistaClear* system
- In-line check valves ("C" above) must **ALWAYS** be used – be sure to install them in the proper flow direction as shown
- The maximum operating pressure for air and water into the *VistaClear* system is 100 psi – preferably 40 psi for water, 40-80 psi for air

"IS-1000-A"  
© 2001 VRG, LLC